MODE D'EMPLOI  
PHYWE AKTIENGESELLSCHAFT  
GÖTTINGEN

GENERATEUR A BANDE  
Alimentation 220 V\*

07642.93

A. BUT

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07642.93

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### A. Purpose

Van de Graaff generators are used to produce high d.c. voltages at comparatively low power. They were developed by the American physicist van de Graaff around 1930 and were the first high voltage sources suitable for research in atomic physics. In the simplified design described here the generator is intended as a demonstration machine to replace the very trouble-prone influence machines. PHYWE van de Graaff generators are safe to work with and do not require any special precautions — a great advantage, particularly when used in schools and colleges.

### B. Working principle

The van de Graaff generator is self-excited. Its working principle can be seen on the accompanying diagram (Fig. 1), which shows only the parts of the generator essential to its functioning. In this illustration, and in those which follow, the various components and parts are referred to by numbers which stay the same throughout.

An endless rubber charging belt (5) runs over two rollers (1 and 2). The lower roller (1) is the driving roller and is made of metal. The excitation roller (2), is made of insulating material (PLEXIGLAS®). Two rows of points (3 and 4) carry electric charges to the belt. The conductor (6), a Faraday cage, collects the charges. The conductor and the excitation roller are carefully insulated from the driving roller. The upper row of points (4) is electri-

cally connected to the conductor. The driving roller (1) is electrically connected to the lower row of points (3) and grounded.

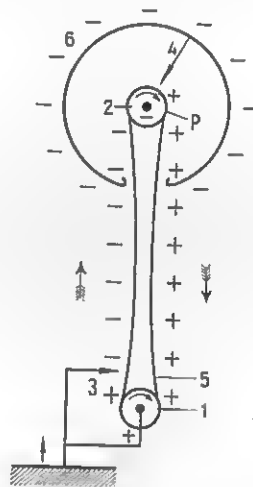
Fig. 1: Principle of van de Graaff generator

The means by which high voltage is generated is as follows: When the driving roller is turned in the arrowed direction, the belt moves and, at point P, breaks contact with the excitation roller. This causes charges of opposite sign to accumulate on the excitation roller and on the belt. Since the machine uses a PLEXIGLAS® roller and a rubber belt, the charge formation is normally negative on the conductor and positive on the belt. It may be found, however, that this charge condition is only obtained after a long period of operation. The positive charge remains on the descending side of the belt which, as it passes the lower row of points, draws negative charges from them by influence effects. These

\* For Cat. Nos. of units designed for operation on other voltages see Section G.

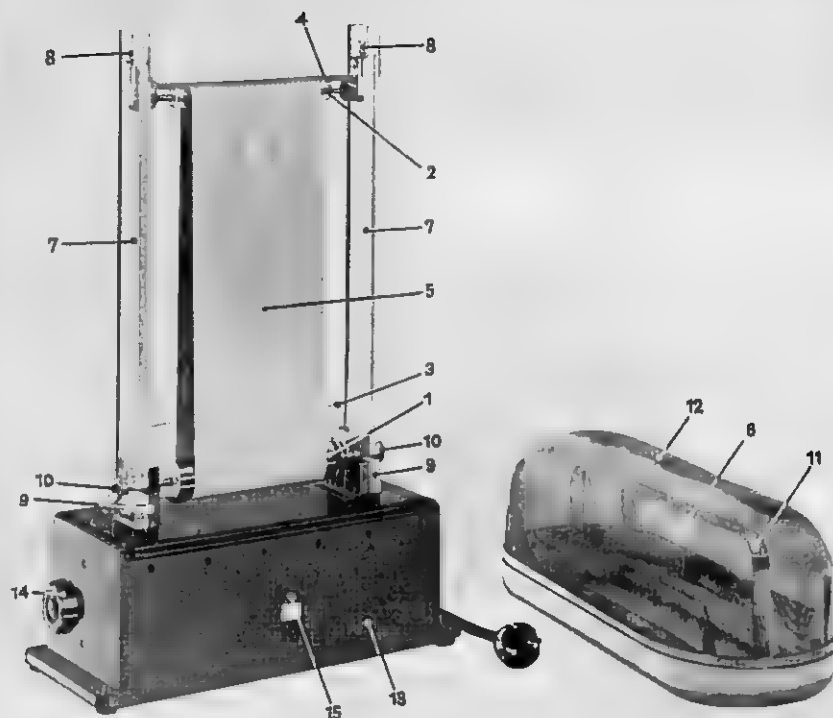
Fig. 2: Details of the generator

charges spray on to the ascending side of the belt, remain there and are carried by the motion of the belt into the conductor. There they cause positive charges to be drawn on to the belt from the upper row of points through influence action. These positive charges are drawn away from the conductor leaving a surplus of negative charges. The process is repeated as the belt continues running. The excess negative charges increase and result in a constant rise in the voltage on the conductor. The voltage obtainable is limited by the insulating properties of the material and by the dielectric strength of the air.

Generally speaking, the polarity of the conductor does not matter for experimental purposes. However, should it be necessary for a particular experiment to have the conductor charged positively, this condition can be obtained by cleaning the generator with solvent (cf. Section E) for a few minutes. The polarity of the conductor can be checked easily with the small neon tube supplied. Hold the tube by its metal cap and bring it up to the conductor. If the charge is negative the electrode nearest the conductor glows. With a positive charge the glow appears on the electrode remote from the conductor.

## C. Generator details

The construction of the generator (Fig. 2) has been kept as simple as possible and closely follows the diagrammatic layout of Fig. 1. The electrically active components are the mirror-polished, chrome-plated driving roller (1), the excitation roller (2) made of PLEXIGLAS®, the lower row of points (3), the upper row of points (4), the belt (5) and the conductor (6).

Rectangular PLEXIGLAS® columns (7) well insulate the excitation roller as well as the upper row of points and the conductor from the driving roller and the lower row of points which is in electrical contact with the roller.

Both rollers run in ball bearings which are pressed on the roller spindles. The ball bearings are fitted into housings which are held in bearing slots. Knurled screws (8) working in grooves in the top ends of the PLEXIGLAS® columns allow the separation between the excitation roller and the driving roller to be adjusted. In this way the mechanical tension of the belt can be varied.

The driving roller is carried in two long slots in the bearing blocks (9). Endwise location in the slots is provided by the knurled nuts (10) at the side.

The high voltage electrode is the conductor (6) which is designed as a Faraday cage. The conductor rests on the heads of the two knurled screws (8) and makes electrical contact via the metal bridge piece (11) fitted inside. The high voltage can be taken off the socket (12) which is fixed to the metal bridge.

The generator is mounted on a moulded plastic base (13) in which is housed the drive unit. This consists of a universal motor with a series resistor operated by the rotary knob (14).

A socket (15) is mounted on each of the long sides of the base. The socket marked with the symbol  $\perp$  is connected to the motor screening, while the socket on the opposite side (not visible in the illustration) marked with the letter S is connected to the bottom row of points, but insulated from the grounded socket. The socket is so insulated to make quite sure that there is absolutely no possibility of connecting a measuring instrument to the line

voltage by accident. The drive is transmitted to the generator through pulleys and a plastic belt ensuring virtually noiseless running.

## D. Placing in service

The generator is packed in a special carton and, except for the conductor, is fully assembled and adjusted. The carton contains the following items:

- 1) van de Graaff generator
- 2) Conductor
- 3) Discharge sphere on stem
- 4) Insulating base for discharge sphere
- 5) Connecting lead, 50 cm long
- 6) Miniature neon tube for polarity check

After removal from the packing material it is advisable to wipe the dust carefully off all parts. While doing this, the excitation roller, driving roller and belt should be treated with special care.

Next place the conductor on the generator so that its metal bridge piece (11) rests on the heads of the knurled screws (8).

After connection to the line voltage the motor can be started and its speed adjusted – and consequently the belt speed – by means of the knob (14) controlling the rheostat. The generator is now ready for operation.

When the stem of the discharge sphere and the insulating base are screwed together, the sphere can then be connected by a lead to the S or  $\perp$  socket. The sphere gives an order-of-magnitude indication of generator voltage on the basis of spark length (1 mm corresponds to about 3 kV). The high voltage is picked off at the socket (12), the opposite pole being the S socket. The  $\perp$  socket must always be grounded. Except when sensitive measuring instruments are being used, the S and  $\perp$  sockets should be connected together by a lead.

## E. Trouble shooting

The van de Graaff generator is so straightforward in construction that mechanical trouble

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P.O.B. 665 · Telephone 61056 · Telex 096808 · Telegrams PHYWE Goettingen

during operation is not likely to occur. However, like all other electrostatic generators, it is sensitive to moisture and may give unsatisfactory results when the atmospheric humidity is high. Any failure occurring under these conditions is in no way attributable to faults in manufacture or design.

When the relative humidity is 70% or higher, satisfactory operation cannot be guaranteed unless special provisions are made. If extreme climatic conditions keep the humidity at a high level constantly, dependable operation can be ensured only if the generator is installed in an air-tight enclosure and the air inside dried with a powerful drying agent such as phosphorus pentoxide ( $P_2O_5$ ).

In the event of the generator's failing to charge up, cleaning will be necessary. For this purpose it must be taken apart.

The actual cleaning process can now be carried out. Using a clean linen cloth soaked in methylated alcohol, first clean the PLEXIGLAS® columns (7). Then clean the driving roller and rub until quite dry. The belt, too, should be cleaned, preferably on both sides, with a cloth soaked in methylated alcohol. The parts cleaned so far should now be replaced in the reverse order to that followed in dismantling. Not until this has been done should the cleaning of the excitation roller proceed as described below, with immediate replacement on completion.

The excitation roller is first wiped off with a linen cloth soaked in methylated alcohol and is then rubbed dry again. A dry cloth is then used to charge up the excitation roller. In other words, the roller is rubbed until it is able to attract small fragments of paper. The

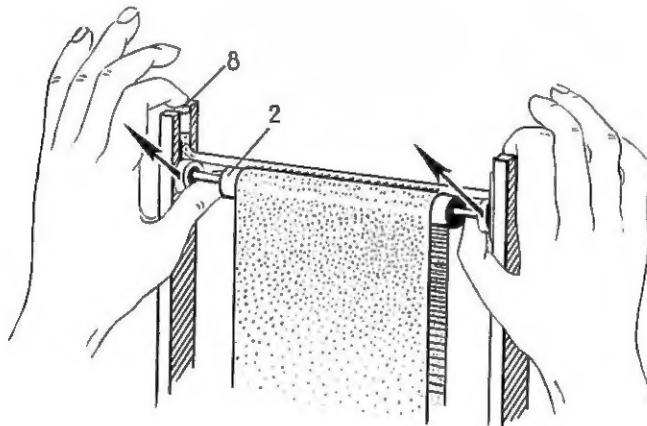
Fig. 3: Lifting out the excitation roller

First, remove the conductor hood and lift the excitation roller out of its bearings at each end simultaneously by drawing upwards at an angle (Fig. 3). At the same time exert an opposing pressure on the knurled screws (8). The excitation roller can now be withdrawn from the belt.

Next, loosen the knurled nuts (10) and pull the driving roller forward out of its bearing slots. The belt can now be taken right out and the plastic driving belt taken off the pulley on the drive shaft. To prevent the driving belt from dropping into the slot in the base of the machine it is advisable to hold it up by passing a pencil or similar object through the loop of the belt.

roller is then passed through the belt and entered into its bearings at either side. When doing this, care should be taken to ensure

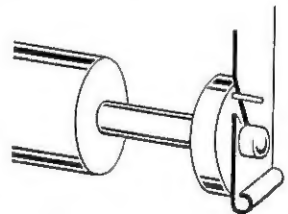
Fig. 4: Entering the roller into its bearing

that the small peg (Fig. 4) on the bearing housing is not introduced into the slot.

In the great majority of cases cleaning the entire machine as described will not be necessary, since the cleaning of the excitation roller alone will produce the desired result. If it is found during operation of the generator that the belt tends to wander off the roller to one side or the other, the necessary correction can be made by the adjustment facility provided by the knurled screws (8) which allow the belt to be adjusted for correct tracking in the middle of the roller.

## F. Performance

The maximum voltage obtainable from the machine depends largely on how clean it is kept and on the climatic conditions. If the generator is looked after properly and if the air is dry, a maximum voltage of 300 kV (a spark length of about 10 cm measured between the discharge sphere and the end of the conductor) can be obtained. At short circuit the current ranges from 10 to 20  $\mu$ A and occasionally higher while the output is 0.3 to 0.5 W.

Apart from keeping it clean, the generator requires little maintenance. After about 100

running hours the ball bearings should be carefully re-oiled with non-gumming oil.

CARE MUST BE TAKEN TO ENSURE THAT NO OIL FINDS ITS WAY ON TO THE BELT, THE EXCITATION ROLLER OR THE PLEXIGLAS® COLUMNS, OTHERWISE THE PERFORMANCE OF THE MACHINE WILL BE GREATLY IMPAIRED.

After it has been in use for an extended time the charging belt should be replaced, since a great deal of ozone is produced resulting in rapid aging of the rubber.

## G. Materials

Designation	Cat. No.	Qty.
van de Graaff generator including discharge sphere on stem with insulating base, connecting lead (50 cm) and miniature neon tube		1
for operation on 110 V a. c.	07642.90	
for operation on 125 V a. c.	07642.91	
for operation on 220 V a. c.	07642.93	
for operation on 250 V a. c.	07642.94	
Belt for van de Graaff generator, replacement	07642.01	1